



- Fits standard 13 x 100 and 16 x 100 collection tubes
- Available with integral filters for prep work
- Color coded


Tip Top™ Dispenser Caps

A safe and cost effective way of dispensing samples in your lab.

The Tip Top is a disposable cone shaped cap that fits all standard 13 or 16mm tubes. A small opening allows you to dispense fluid without removing the Tip Top cap so you can prepare smears on slides or dispense into reagent well or micro centrifuge tubes.

- **Convenient** - just tilt the tube and squeeze to dispense one drop
- **Safe** - you don't need to open the tube and pipet out of the tube so spillage is reduced along with potential exposure to infectious disease
- **Low Evaporation** - A flexible cap is available to cover the Tip Top so you can store samples in the refrigerator.
- **Color Coded** - available in three colors for easily lab organization

The Tip Top cap can be used with almost any liquid. Made from chemically resistant polyethylene they can dispense serum, reagents, plasma, whole blood and CSF. A filtered version is available for use as a serum separator or reagent blender.

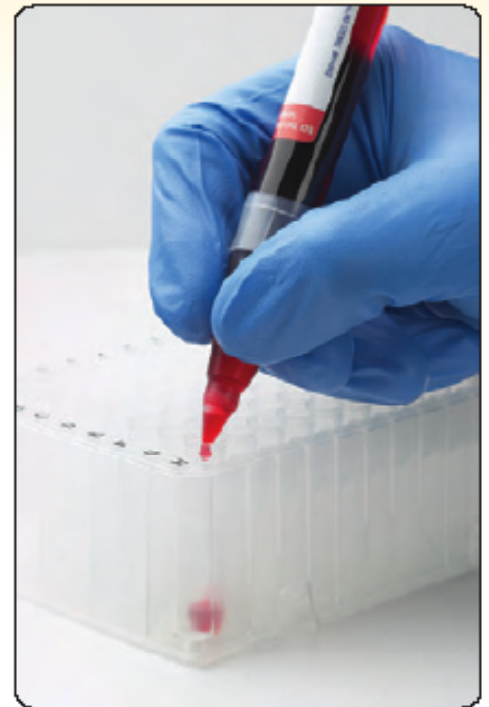

Cat No.	Description	Unit	Case
<i>Tip Tops for 13mm tubes</i>			
5501-800	Natural, 100/resealable bag	1000	5000
5503-804	Red, 100/resealable bag	1000	5000
5505-807	Blue, 100/resealable bag	1000	5000
<i>Tip Tops for 16mm tubes</i>			
5502-800	Natural, 100/resealable bag	1000	5000
5504-804	Red, 100/resealable bag	1000	5000
5506-807	Blue, 100/resealable bag	1000	5000
<i>Tip Tops with Integral Filters</i>			
5511-800	Natural, Fits 13mm tubes, 1000/bag	1000	5000
5512-800	Natural, Fits 16mm tubes, 1000/bag	1000	5000
<i>Closures for Sealing Tip Top Caps</i>			
5515-804	Closure, Red, 1000/bag	1000	5000